

RENAISSANCE[®]



Nice Care
Health Products Limited



www.nicecare.uk

MAGNETIC THERAPY DEVICES CATALOGUE

- ✚ for home and professional use
- ✚ assists you in selecting the most suitable device to suit your needs
- ✚ includes a map with application points on your body



TABLE OF CONTENTS

Magnetic therapy as a healing procedure	4
Renaissance brand introduction	6
How to select the most suitable magnetic therapy device to suit your needs	10
DUO FORTE generator	20
MANIC applicator	22
PIANO applicator	25
LETTINO applicator	28
CASSO applicator	31
CAVO applicator	34
VIERA applicator	39
PRIMO C set	42





PRIMO M set	44
ANELLI set	46
ABOLITO set	48
MELIO set	50
MEDICO set	52
LETTINO set	54
VIVA set	56
VIVA PLUS set	58
SIECCO set	61
VIERA set	63
When magnetic therapy should not be used or caution exercised	65
Application points to be used during magnetic therapy	66
Contact details	69






MAGNETIC THERAPY AS A HEALING PROCEDURE

Scientific research and our experience with Renaissance devices indicate pulsed magnetic therapy has a number of healing properties. These can be summarized as follows: Magnetism induced by pulsed low-frequency magnetic field is an efficient alternative to

alleviates pain	analgesic effect
removes muscle tension	myorelaxant effect
makes blood vessels wider and stabilizes blood pressure	vasodilatation effect
treats swelling	anti-inflammatory and anti-rheumatic effect
stimulates growth and construction of bone tissues	stimulating effect
eliminates harmful substances and accelerates metabolism	detoxicating effect
regenerates cells and improves blood supply	regenerating effect
accelerates healing processes	healing effect
boosts immunity	immunity effect

Magnetism induced by pulsed low-frequency magnetic field **is an efficient alternative to currently overused pharmacological treatment**. One of the greatest assets of this type of treatment is that it is natural. As opposed to pharmacological treatment, magnetism **does not create addiction**.

In addition to the effects given above, magnetic therapy is an incredibly efficient way to treat **swelling** (oedema) at home. A magnetic field also makes small blood vessels wider as well as improves blood supply, properties that are used for treating **diabetes**. Additionally, magnetic therapy is highly efficient in treating **migraines** and **ringing in the ear** (tinnitus) and is also helpful in treating a number of psychological disorders (such as sleep disorders) and autoimmune diseases (**Bechterew's disease, rheumatism**). Renaissance provides clients



with specially designed magnetic therapy sets that can also be used during **post-operative convalescence**, to accelerate **healing of wounds** as well as treat **paradontosis, arthritis and arthrosis, sciatica** and other **musculoskeletal disorders**. The sets have also proved to be a highly efficient method of dealing with **lifestyle diseases**.

Magnetic therapy is not to be used in older and unhealthy organisms only. Sportspeople know Renaissance magnetic therapy products as they frequently use them for the purposes of physical medicine and rehabilitation as well as to relieve sore muscles after strenuous exercise.

Pulsed magnetic therapy has become a very popular alternative treatment method in the course of the last decade. Experts currently consider magnetic field treatment (magnetic therapy) **to be a safe and highly efficient physical treatment method**. It is often referred to as “the treatment method of the 21st century”. It should be noted that what is involved is not the use of simple static magnets (such as wristbands, bandages or necklaces) but highly efficient low-frequency pulsed magnetic field.

The healing and regenerating properties of electromagnetic energy had until recently only been used in facilities such as hospitals and professional physical medicine centres. Many years of experience and social changes and development, however, have gradually transformed pulsed magnetic therapy into a treatment method available for home use as well.



Pulsed magnetic therapy has successfully gained ground in medical fields such as **orthopedy, sports medicine, balneology, physical therapy and rehabilitation, pediatry, geriatrics, dermatology, neurology, stomatology, urology, internal medicine** and **veterinary medicine**.



RENAISSANCE[®]

THE RENAISSANCE BRAND

Long and strenuous research had to be carried out before Renaissance magnetic therapy sets entered production. The result is magnetic therapy that fulfils the requirements of both **professional and home use**. It provides efficiency and is very simple to use during **home application** and for the purposes of **self-therapy**.

One of the assets of the Renaissance system is it was originally intended to be used in hospitals and in clinical practice. This means it is designed to be highly durable, boast high quality materials, and withstand repeated, intensive use. Home users appreciate the system is incredibly user-friendly, safe and offers considerable versatility in terms of what conditions it can treat.

The Renaissance brand uses its own generator combined with a variety of applicators to suit each client's specific needs. The Renaissance magnetic pulse generator is optimized to create the most suitable frequency range, making it possible for users to regulate both the frequency and the length of application. In contrast, applicators are able to “transport” the healing magnetic field to where it is most needed. Renaissance applicators can be used locally to treat specific points as well as larger body sections during bedside applications. And our range also includes “depth” circular applicators.

Renaissance magnetic therapy products are quickly expanding throughout Europe and overseas. What also demonstrates their success is that existing clients purchase new Renaissance applicators to substitute or supplement the ones they purchased earlier.

The most recent Renaissance healing sets for low-frequency pulsed magnetic therapy are **modern certified medical equipment** designed for home use as well as for professional physiotherapists. They use state-of-the-art technologies, are highly efficient and user-friendly, and boast an ergonomic design. As a result, they fulfil all the challenging requirements contemporary medicine might have.

MAJOR REASONS TO CHOOSE RENAISSANCE MAGNETIC THERAPY PRODUCTS



Extended warranty

We stand for the quality of our devices. That's why main control units of our products come with **an extended warranty**. **Renaissance magnetic therapy apparatuses are compliant with all the norms applicable to medical equipment.** Their parameters and function are regularly checked.



Life-long care guarantee

We offer all our clients **professional maintenance services throughout the life of the machine**. Additionally, we are happy to **advise you free of charge** whenever necessary as well as refer you to a medical doctor to answer your questions.



We represent quality

The best way to demonstrate quality of medical equipment is to show it really works. We have received **a large number of recommendations from our satisfied customers**. Our magnetic therapy devices boast exceptional quality and reliability. Renaissance products' high and **certified quality** has granted the brand a key position on the Czech, Slovak as well as other markets in and outside Europe.



Tradition since 1995: production, research, development

Our establishment, which is found in the industrial park of Brno-jih and which covers an area of over 34 acres, **includes our own production, administrative and storage spaces.**

We have been carrying out research and developing and producing complete Renaissance devices since 1995. Our company is able to address the needs of our clients both in terms of the number of devices we can produce and by offering complex advice and consulting services related to use and application of magnetic therapy. You can rely on us.



Made in the Czech Republic

Renaissance products are, in their entirety, manufactured in the Czech Republic to be subsequently distributed to the whole world. **Over twenty thousand satisfied customers from the Czech Republic recommend Renaissance magnetic therapy products.** The manufacturer has been a proud member of the Association of Manufacturers and Suppliers of Medical Devices existing under the aegis of the Ministry of Health of the Czech Republic.



Hand-made

Applicators used with Renaissance magnetic therapy devices are **hand-sewn** in our own sewing workroom. They are designed to withstand frequent and long-term abuse. Their surface is easy to wash. And their production quality has been acknowledged by users and manufacturing process professionals.



For both home use and professionals

Renaissance magnetic therapy devices have been used, on a daily and long-term basis, by **professionals** in a number of hospitals, private healthcare facilities, physical medicine and rehabilitation centres and other specialized facilities. They are regularly used by doctors, physiotherapists, and masseurs as well as in wellness and fitness centres. Many of them are currently being purchased to be **used at home**, where they help individuals and whole families to deal with a number of health issues as well as during convalescence and recovery. Magnetic therapy devices are also often used for preventative reasons. In addition, Renaissance devices are sought after by **sportspeople** and may also be used to efficiently treat **animals**.



A product for the whole family

When used at home, Renaissance magnetic therapy devices offer **a wide variety of applications**. They can be used to treat a number of health issues, improve mobility, alleviate pain, improve physical condition, boost immunity system, and accelerate healing processes and make them more efficient. Additionally, the devices may be used for preventative purposes as well as to help you to relax. Magnetic therapy can be enjoyed by **children, adolescents, adults, and elderly people**. And your **pet** can experience its positive effects too.




User-friendly

Everyone is able to connect and control Renaissance magnetic therapy devices. Each of these devices comes with a comprehensive user's manual and documentation. We are happy to advise you and assist you with whatever needs you might have. And, if necessary, we will personally demonstrate to you how the device works.

SELECTING THE MOST SUITABLE MAGNETIC THERAPY DEVICE TO SUIT YOUR SPECIFIC NEED

Each magnetic therapy device consists of a pulse generator and an applicator (or a combination of applicators).

A pulse generator serves as a control unit where required treatment program can be selected. In contrast, **an applicator**, which is connected with the generator, is to be placed on a specific part of your body which you want magnetic therapy to treat. Once application has finished, the device automatically produces a sound signal. Renaissance magnetic therapy devices can be controlled in a very easy and user-friendly way.

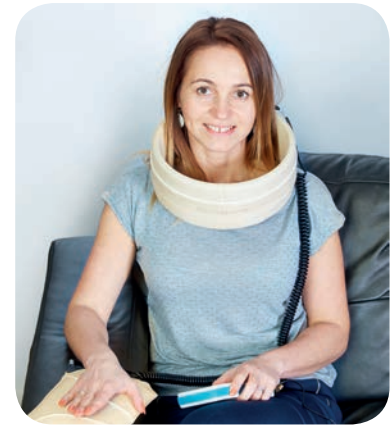

Applicator



Generator

A list of applicators we offer is provided on pages 22-39 of this catalogue. **Please, select applicators based on what painful and affected body parts you need to treat, not based on recommendations for specific illnesses and health issues.**



The following table might help you to select the right kind of applicator:

Diagnosis	Recommended applicator	Where to apply
Bartholin's abscesses	small, flat (MANIC)	lower abdomen
Cold caused by an allergy	small, circular (CAVO)	head, cervical spine
Alzheimer's disease	small, circular (CAVO) and large, flat (PIANO, LETTINO, CASSO)	head, cervical spine
Amputation stumps	small, circular (CAVO)	affected area
Angina pectoris, IHD	small, circular (CAVO) and flat (MANIC, PIANO, LETTINO, CASSO)	chest from front and back
Tonsillitis, nasopharyngeal infections	small, circular (CAVO)	neck
Arthritis in patients suffering from psoriasis	small, circular (CAVO)	affected area
Arthrosis of joints	small, circular (CAVO)	affected area
Knee arthrosis	small, circular (CAVO)	knee area
Hip arthrosis	small, circular (CAVO) or large, circular (VIERA)	hip area
Arrhythmia	small, flat (MANIC)	chest area
Aseptic necrosis of bones	small, circular (CAVO)	affected area
Asthma (bronchial asthma)	small, circular (CAVO) or small, flat (MANIC)	thymus area
Atopic dermatitis	small, flat (MANIC)	affected place + that spinal area which is the source of innervation
Optic nerve atrophy	small, circular (CAVO)	head area
Bechterew's disease	large, flat (PIANO, LETTINO, CASSO)	spinal area
Varicose ulcers	flat (MANIC, PIANO, LETTINO, CASSO) or small, circular (CAVO)	affected area




Diagnosis	Recommended applicator	Where to apply
Back and neck pain	large, flat (PIANO, LETTINO, CASSO)	spinal area
Feeling of pins and needles in limbs (paresthesia)	large, flat (PIANO, LETTINO, CASSO) and circular (CAVO, VIERA)	affected area and cervical or lumbar spine
Bronchitis	small, circular (CAVO) or small, flat (MANIC)	thymus area
Cellulitis	large, flat (PIANO, LETTINO, CASSO)	affected area
Spinal chord degeneration	small, flat (MANIC)	spinal chord area
Degenerative retinoschisis	small, circular (CAVO)	head area
Detoxication program	large, flat (PIANO, LETTINO, CASSO), circular (CAVO, VIERA)	whole body
Polio in children	small, circular (CAVO) and large, flat (PIANO, LETTINO, CASSO)	head and spinal area
Diabetes	small, flat (MANIC)	pancreas and liver area
Diabetic leg – promoting blood circulation, pins and needles	large, flat (PIANO, LETTINO, CASSO)	affected area
Gout	small, circular (CAVO)	affected area
Dupuytren's contracture	small, flat (MANIC)	palm area
Endometriosis	large, flat (PIANO, LETTINO, CASSO) and small, circular (MANIC)	lumbar and lower abdominal area
Endoprostheses	small, circular (CAVO) or large, circular (VIERA)	joint area
Erectile dysfunction	small, flat (MANIC) or circular (CAVO, VIERA)	genital area



Diagnosis	Recommended applicator	Where to apply
Cervical erosion	small, flat (MANIC)	lower abdominal area
Glaucoma	small, circular (CAVO)	head area
Gynaecological surgeries	small, flat (MANIC)	lower abdominal area
Haemorrhoids	small, flat (MANIC)	rectal area
Larynx, vocal chords, hoarse voice	small, circular (CAVO)	neck
Hypothyreosis	small, circular (CAVO)	laryngeal area
Immunity system boost	small, circular (CAVO)	thymus area
Sciatica	large, flat (PIANO, LETTINO, CASSO)	spinal area
Liver damage caused by hepatitis or exposure to toxins	small, flat (MANIC)	liver area
Carpal tunnel syndrome	small, flat (MANIC) and circular (CAVO, VIERA)	wrist area and cervical spine
Postmenopausal syndrome	small, flat (MANIC)	lower abdominal area
Tick-borne encephalitis	small, circular (CAVO) and large, flat (PIANO, LETTINO, CASSO)	head and spinal area
Breast-feeding and lactation suppression	large, flat (PIANO, LETTINO, CASSO) or small, flat (MANIC)	chest area
Blocked cervical spine	small, circular (CAVO)	cervical spine
Varicose veins	large, flat (PIANO, LETTINO, CASSO)	head area
Minimal brain dysfunction in children (MBD)	small, circular (CAVO)	oblast hlavy
Psoriasis	small, flat (MANIC) and circular (CAVO, VIERA)	affected area



Diagnosis	Recommended applicator	Where to apply
Lyme disease	large, flat (PIANO, LETTINO, CASSO) and small circular (CAVO)	spinal chord and head area
Massage	large, flat (PIANO, LETTINO, CASSO)	whole body
Ménière's disease	small, circular (CAVO)	head area
Painful, irregular menstruation	small, flat (MANIC)	lower abdominal area
Migraine	small, circular (CAVO)	head area
Mononucleosis	small, flat (MANIC) and circular (CAVO, VIERA)	liver and cervical spine area
Brain stroke – consequences	small, circular (CAVO)	head area
Myopathy	large, flat (PIANO, LETTINO, CASSO)	spinal area
Chronic obstructive pulmonary disease (often caused by smoking)	large, flat (VIERA) or flat (MANIC, PIANO, LETTINO, CASSO)	chest area
Infertility (reproductive system disorders)	large, flat (VEIRA)	loins and lower abdominal area
Insomnia (sleeping disorders)	small, circular (CAVO)	head area
Odontogenic osteomyelitis	large, flat (VIERA) or small, circular (CAVO)	affected area
Lymphatic vessel diseases	large, flat (PIANO, LETTINO, CASSO)	whole body
Artery diseases	large, flat (PIANO, LETTINO, CASSO)	whole body



Diagnosis	Recommended applicator	Where to apply
Vein diseases	large, flat (PIANO, LETTINO, CASSO) or small, flat (MANIC)	affected area
Osteoporosis	large, flat (PIANO, LETTINO, CASSO) or circular (CAVO, VIERA)	affected area
Spurs	small, flat (MANIC)	foot and heel area
Swelling	flat (MANIC, PIANO, LETTINO, CASSO) or circular (CAVO, VIERA)	affected area
Heartburn, reflux disease	small, flat (MANIC) or circular (CAVO, VIERA)	stomach and thoracic spine
Paradontosis	small, flat (MANIC) or circular (CAVO, VIERA)	head area
Parkinson's disease, stiffness	large, flat (PIANO, LETTINO, CASSO) or small, circular (CAVO)	spinal and head area
Herpes zoster	large, flat (PIANO, LETTINO, CASSO)	affected area
Post-operative spine convalescence	large, flat (PIANO, LETTINO, CASSO)	spinal area
Spinal apophyses	large, flat (PIANO, LETTINO, CASSO)	spinal area
Unspecified pulmonary diseases	large, circular (VIERA)	chest area
Postpartum uterine regression	large, circular (VIERA)	loins and lower abdominal area
Urinary incontinence in children	small, flat (MANIC)	lower abdominal area



Diagnosis	Recommended applicator	Where to apply
Post-operative conditions (healing)	large, flat (PIANO, LETTINO, CASSO) or circular (CAVO, VIERA)	affected area
Burns	small, flat (MANIC) or large, flat (PIANO, LETTINO, CASSO)	affected area
Limb injuries	small, circular (CAVO)	affected area
Peripheral neuropathy	small, circular (CAVO)	head area
Prostate – inflammation, enlargement	large, circular (VIERA) or flat (MANIC, PIANO, LETTINO, CASSO)	lower abdominal area
Nipple rhagades	small, flat (MANIC) or large, flat (PIANO, LETTINO, CASSO)	nipple area
Rheumatoid arthritis	large, flat (PIANO, LETTINO, CASSO)	spinal area
Cornea – injuries, operations, eye infections	small, circular (CAVO)	head area
Multiple sclerosis	large, flat (PIANO, LETTINO, CASSO) and small, circular (CAVO)	head and spinal area
Sinusitis – sinus infections	small, circular (CAVO)	head area
Scleroderma – neurological complications	small, circular (CAVO)	head area
Tennis elbow	small, circular (CAVO) or flat (MANIC, PIANO, LETTINO, CASSO)	elbow area
Tinnitus	small, circular (CAVO)	head area



Diagnosis	Recommended applicator	Where to apply
Pulmonary tuberculosis (chronic)	large, circular (VIERA)	chest area
Ulcerative colitis	large, circular (VIERA)	abdominal, loin and hip areas
Spinal chord injuries	small, flat (MANIC)	spinal chord area
Low back pain (lumbago)	large, flat (PIANO, LETTINO, CASSO)	spinal area
Relaxation of muscles (myorelaxation)	large, flat (PIANO, LETTINO, CASSO)	affected area
Herniated spinal disc	large, flat (PIANO, LETTINO, CASSO)	spinal area
Hypertension (high blood pressure)	large, flat (PIANO, LETTINO, CASSO)	chest area
Constipation	large, flat (PIANO, LETTINO, CASSO) and small, flat (MANIC)	abdominal area
Uterine and endometrial inflammation	small, flat (MANIC)	lower abdominal area
Bone marrow inflammation (osteomyelitis)	small, circular (CAVO)	affected area
Kidney inflammation, colic	large, circular (VIERA)	loin area
Spinal root inflammation (polyradiculoneuritis)	small, flat (MANIC)	spinal cord area
Urinary tract inflammation	small, flat (MANIC)	lower abdominal area
Vaginal inflammation and discharge	small, flat (MANIC)	lower abdominal area
Inflammation in pancreas (acute)	large, circular (VIERA)	transition between thoracic and lumbar spine



Diagnosis	Recommended applicator	Where to apply
Inflammation in pancreas (chronic)	large, circular (VIERA)	transition between thoracic and lumbar spine
Inflammatory disease of the middle ear (otitis media)	small, circular (CAVO)	head area
Ulcerative colitis, Crohn's disease	large, circular (VIERA) or flat (MANIC, PIANO, LETTINO, CASSO)	abdominal area
Ovarian inflammation and inflammation in fallopian tubes	small, flat (MANIC)	lower abdominal area
Fractures	small, circular (CAVO) or flat (MANIC, PIANO, LETTINO, CASSO)	affected area
Fractures – complicated healing	small, circular (CAVO) or large, circular (VIERA)	affected area
Frozen shoulders	small, circular (CAVO)	shoulder area
Gastric ulcers	large, circular (VIERA)	stomach area

HOW TO USE APPLICATORS – RECOMMENDATIONS



While using an applicator, please, bear in mind **a spine always works as a whole**. That means if you're having problems in the lower back, it is still necessary to treat the whole spine.

If you're experiencing pain a magnetic therapy device must be applied to the sore body area as well as to that spinal area which constitutes the source of the sore area's innervation


(innervation = the nerve supply, usually to a specific part of the body, during which a specific nerve is used). It is therefore recommended also to apply a magnetic device to your cervical spine although what you want to treat is your arms or hands. Similarly, if you need to treat your legs, it is highly advisable to apply the device to your lower back at the same time.

The best results in the shortest time possible are always achieved by whole-body application.

A circular applicator is the best choice for joint and cervical spine issues as its direct-current magnetic field is able to penetrate the whole depth of the tissue.

Please, consult the pages containing detailed information about the specific applicators. This is where you will also find a number of practical demonstrations of how the applicators are designed to be used. An application point map is provided on page 67. This map gives you a better overview of how the specific applicators can be used.

Each applicator description also contains pictograms (graphic signs) that help you to choose the most suitable applicator to address your needs. These pictograms represent the body part the specific applicator can be applied to.

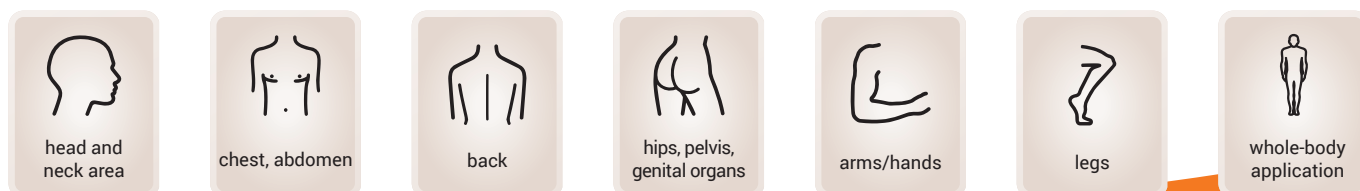

DUO FORTE GENERATOR



DUO FORTE generator description and use:

DUO FORTE is the most up-to-date type of generator with dozens of special programs. Each disease has its own special program combining a number of frequencies and pulse shapes. For example, arthrosis is a disease that affects not only the articular cartilage but also the bone underneath (bone changes, apophyses) and the articular capsule. That is why the specific program has been developed so as to treat the cells in each of these three kinds of tissue, ensuring better therapy results and faster pain relief.

A user of this generator can set intensity and application length as well as choose a different program for each applicator. Frequently-used programs can be stored in the “Favourites” folder to be subsequently restored without the user having to set the desired application length or magnetic field strength again.

DUO FORTE generator parameters:

- a few dozens of special programs
- analgesic programs
- regenerating programs
- two independent output channels that make it possible to run two different procedures at the same time
- regulated optimum frequency (1-100 Hz range)
- easy to use
- ergonomic design








MANIC APPLICATOR



MANIC applicator parameters:

- flat shape
- weight of 455 grams
- dimensions: 170 x 210 mm

Body areas MANIC applicator is the most suitable for:

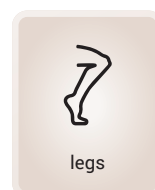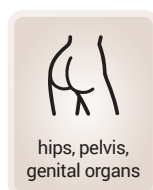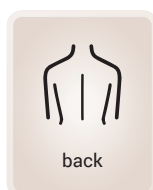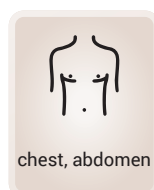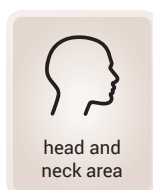


MANIC applicator description and use:


This is a small applicator designed for targeted application to areas such as: ears, teeth, cervical spine, smaller joints, regeneration of inner organs (liver, stomach, uterus, prostate, etc.).

Further information on how MANIC applicator can be used is available in the PRIMO M set description on page 44.

MANIC applicator demonstration:

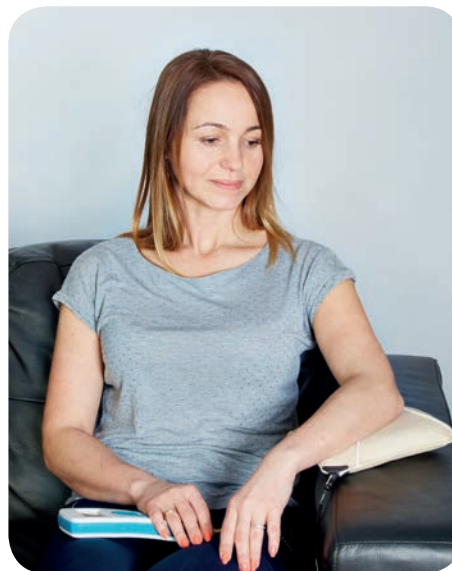

MANIC applicator demonstration:

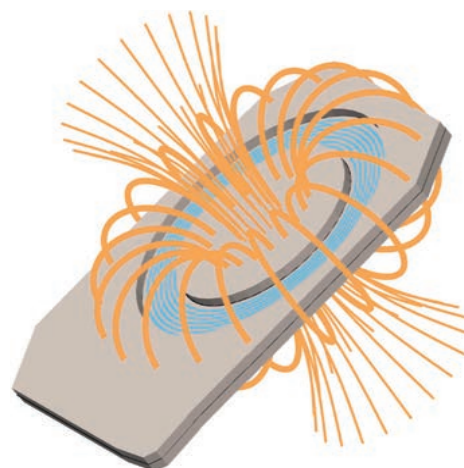

*MANIC applicator magnetic
field lines of force*





PIANO APPLICATOR



PIANO applicator parameters:

- lat shape
- weight of 860 grams
- dimensions: 460 x 490 mm

Body areas PIANO applicator is the most suitable for:

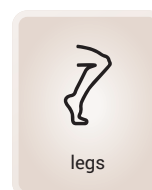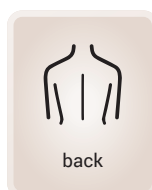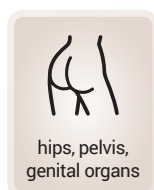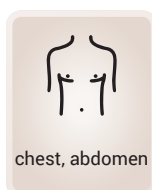






PIANO applicator description and use:

This is an applicator primarily designed to treat the back. However, it's also perfectly suitable to treat organs in the thoracic (chest) cavity and the abdominal cavity. Can also be used for limbs and large joints.


PIANO applicator demonstration:

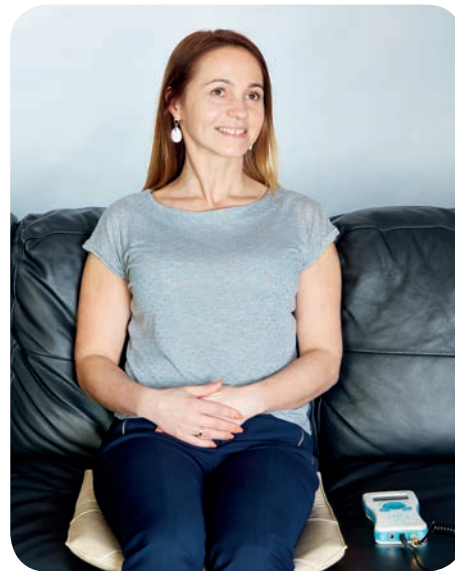

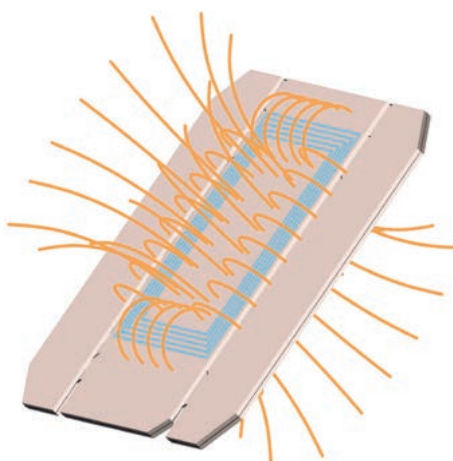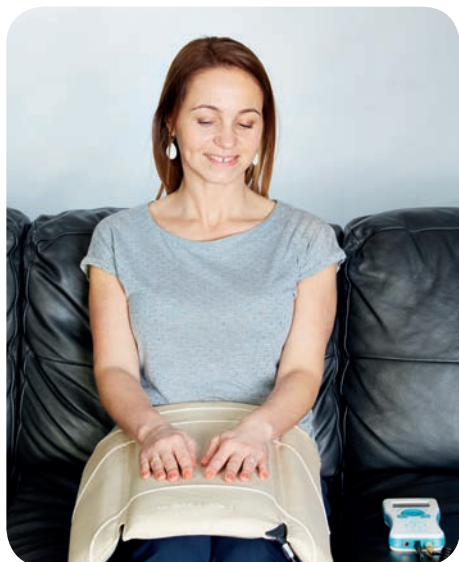

*PIANO applicator magnetic
field lines of force*





LETTINO APPLICATOR

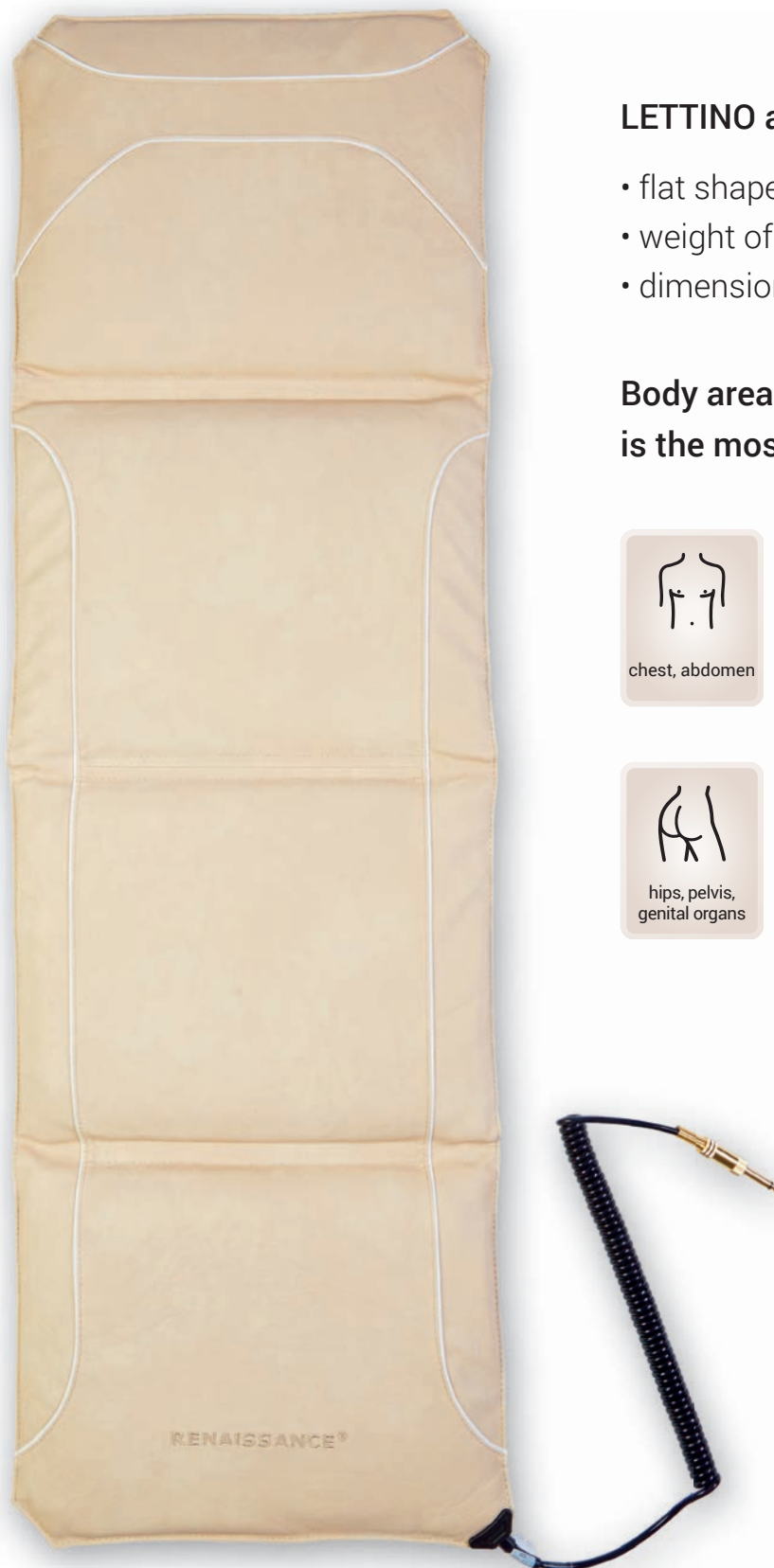

LETTINO applicator parameters:

- flat shape, consists of 4 pieces
- weight of 2,275 grams
- dimensions: 310 x 1,010 mm

Body areas LETTINO applicator is the most suitable for:

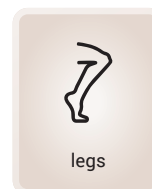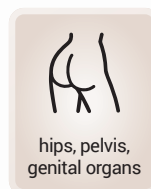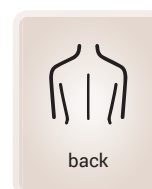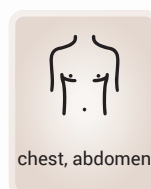


LETTINO applicator description and use:

This is a medium-sized “bed-type” applicator. It is suitable to be applied to the whole spinal area including hips and head. It is also an incredibly efficient therapy device for lower limbs and their joints and is therefore highly recommended to people suffering from diabetes.

LETTINO applicator demonstration:




LETTINO applicator demonstration:

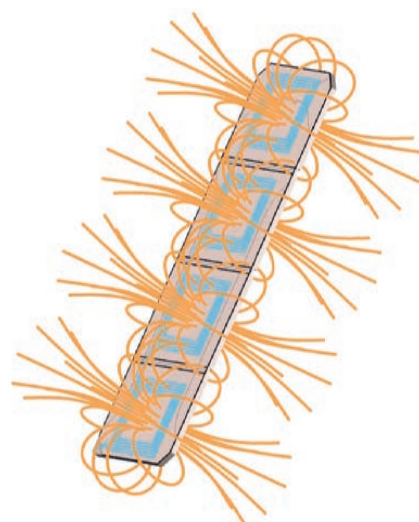

*LETTINO applicator magnetic
field lines of force*


CASSO APPLICATOR

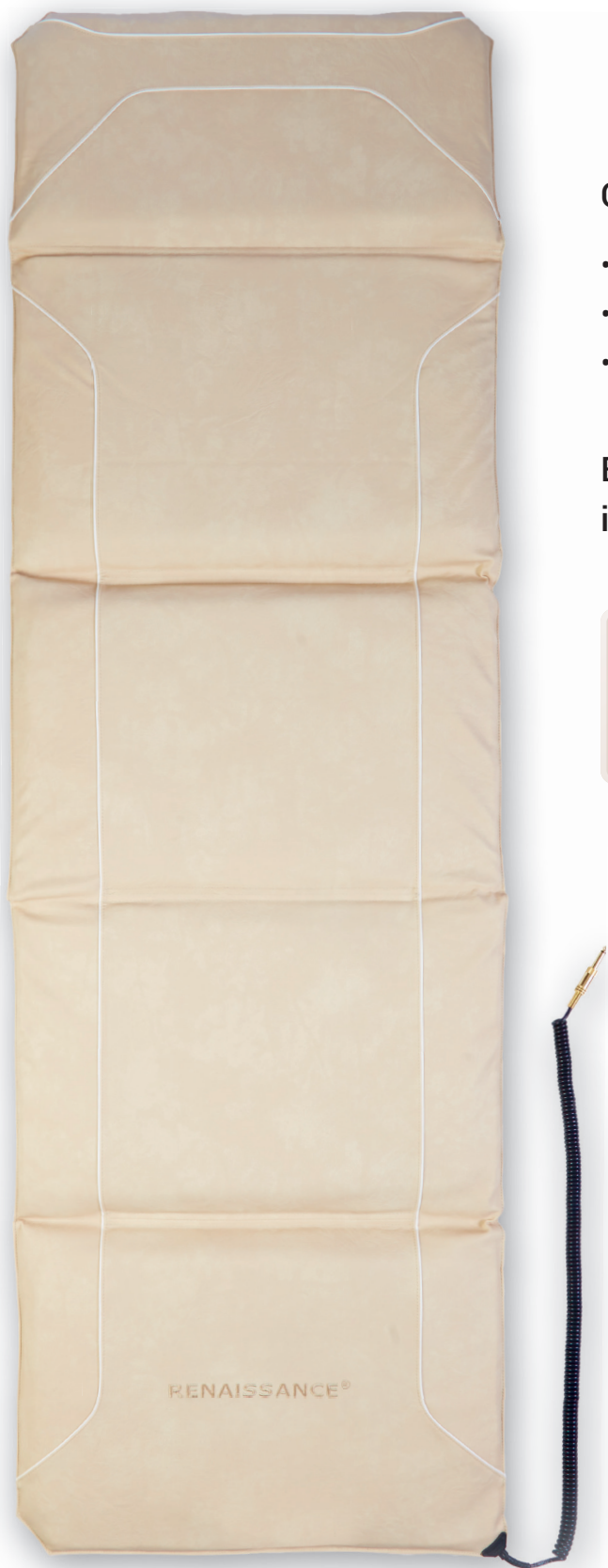

CASSO applicator parameters:

- flat shape, consists of 5 pieces
- weight of 5,980 grams
- dimensions: 550 x 1,730 mm

**Body areas CASSO applicator
is the most suitable for:**

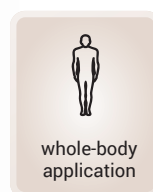


CASSO applicator description and use:

This is the largest “bed-type” applicator. It is suitable for whole-body regeneration.

CASSO applicator demonstration:



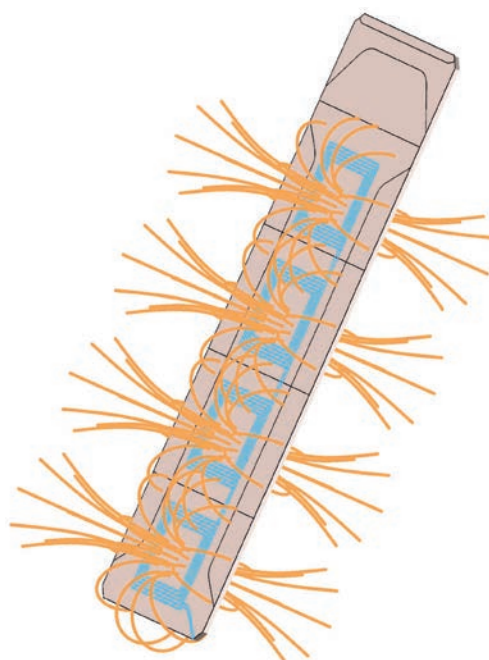

*CASSO applicator magnetic
field lines of force*


CAVO APPLICATOR



CAVO applicator parameters:

- circular shape
- weight of 1,330 grams
- diameter of 260 mm

Body areas CAVO applicator is the most suitable for:

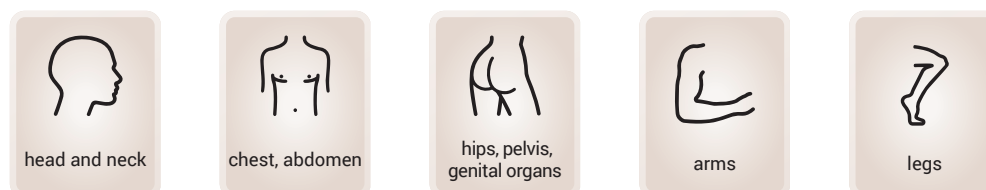

CAVO applicator description and use:

This is the most suitable applicator for people suffering from diseases affecting joints, cervical spine, headache, toothache, ringing in the ear, some neurological disorders, etc. A homogeneous magnetic field in a hollow coil ensures the entire tissue depth at application point regenerates.

Further information on how CAVO applicator can be used is available in the PRIMO C set description on page 44.

CAVO applicator demonstration:






CAVO applicator demonstration:



CAVO applicator demonstration:



CAVO applicator demonstration:

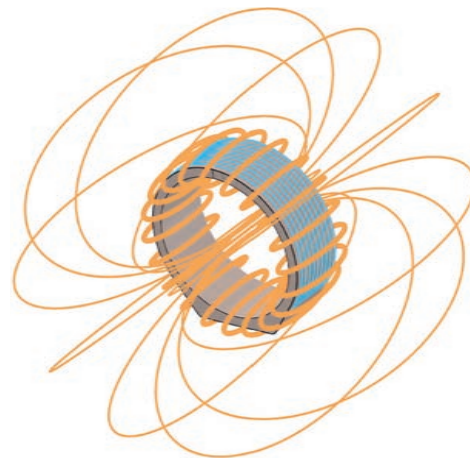

*CAVO applicator magnetic
field lines of force*







VIERA APPLICATOR



VIERA applicator parameters:

- circular shape
- weight of 9.2 kilograms
- diameter: 650 mm

Body areas VIERA applicator is the most suitable for:

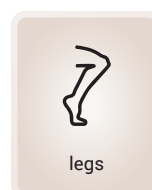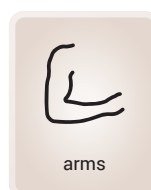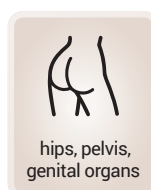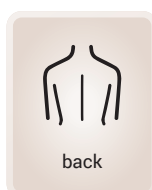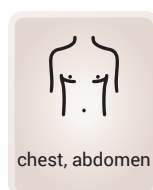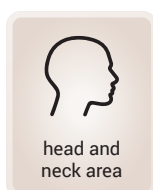




VIERA applicator description and use:

This is a large circular applicator with a diameter of 65 centimetres. It is primarily suitable for professional spa and physical medicine/rehabilitation facilities as well as physiotherapists' offices and, last but not the least, healthcare providing homes and day-care centres for elderly and disabled people. This applicator makes it possible for a patient to be placed inside a homogeneous (same-current) pulsed magnetic field that evenly penetrates the entire depth of the tissues being treated. The applicator's regenerating, analgesic and anti-inflammatory effects are fully indicated in patients suffering from arthrosis in hip and knee joints, rheumatoid arthritis, backache regardless of origin, osteoporosis and a number of other diseases. While pulsed magnetic therapy is an efficient treatment method to fight all the foregoing diseases, it is especially recommended to treat osteoporosis, a disease whose riskiest symptom is decreased bone density.

VIERA applicator demonstration:



VIERA applicator demonstration:

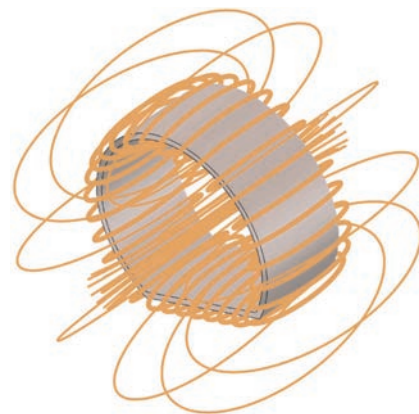

*VIERA applicator magnetic
field lines of force*


PRIMO C SET



PRIMO C set includes:

- a circular CAVO applicator with a diameter of 260 mm
- a two-output pulse generator (DUO FORTE)
- power supply
- magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas PRIMO C set is the most suitable for:

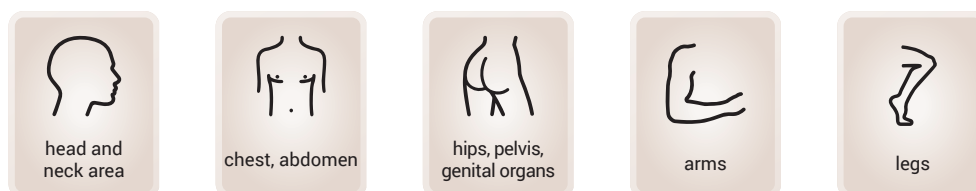

PRIMO C set description and use:

CAVO, an applicator included in this set, has a circular shape with a diameter of 26 cm. It generates homogenous magnetic field with an effective range of approximately 30-40 cm that spreads from the applicator to all sides along a longitudinal axis crossing its centre. This means if we use this applicator to treat pain and dysfunctions in cervical spine, we also regenerate the thoracic (chest) part of spine, respiratory tract and thymus (immunity) as well as the entire head area (migraine, frequent pain), oral cavity (toothache, paradontosis, inflammations), upper respiratory tract, paranasal cavities, eyes (glaucoma, conjunctivitis) and ears (ringing in the ears, middle ear inflammations).

PRIMO C set, which includes a CAVO applicator, is also an efficient alternative to treat a number of dysfunctions in musculoskeletal system, especially joints. As was described above, a CAVO applicator generates a homogenous magnetic field that evenly penetrates the entire tissue depth. As a result, its regenerating, anti-inflammatory and analgesic effect is experienced in all parts of the joint being treated – cartilage, ligament and bone.

PRIMO C set applicator



A circular applicator CAVO

- diameter: 260 mm

Further information about CAVO applicator is available on page 34







PRIMO M SET



PRIMO M set includes:

- a flat **applicator MANIC**, dimensions: 170 x 210 mm
- a **two-output pulse generator (DUO FORTE)**
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas PRIMO M set is the most suitable for:

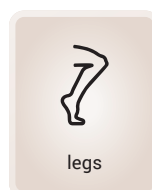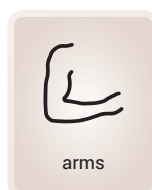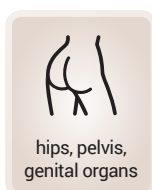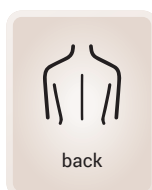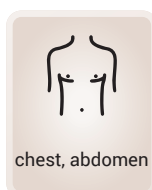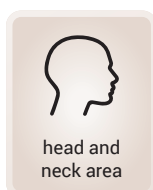

PRIMO M set description and use:

PRIMO M set contains a single small flat MANIC applicator. Despite its humble dimensions this treatment set can help with a surprisingly wide range of health issues.

The MANIC is a small applicator designed to be used for intensive highly targeted therapy and regeneration of various organs and tissues. It efficiently relieves pain in cervical spine as well as gets patients rid of muscular spasms in the area. It is of great help in patients suffering from ringing in the ear (tinnitus, Ménière's disease) and otitis media. And it is also highly recommended for patients with toothache and paradontosis. The size of MANIC applicator makes it possible for it to generate an intensive magnetic field which can alleviate a number of kinds of pain. Its analgesic and regenerating effect is highly appreciated by people suffering from arthrosis or rheumatoid arthritis in knee and hip joints as well as patients with tennis elbows and carpal tunnel syndrome.

PRIMO M set also efficiently treats haemorrhoids and prostate issues and is excellent for targeted regeneration of liver, stomach (gastric ulcers), and other internal organs.

PRIMO M set applicator:



A flat MANIC applicator

- dimensions: 170 x 210 mm

Further information about MANIC applicator is available on page 22


ANELLI SET



ANELLI set includes:

- a circular **applicator CAVO**, diameter: 260 mm
- a flat **applicator MANIC**, dimensions: 170 x 210 mm
- a **two-output pulse generator (DUO FORTE)**
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas ANELLI set is the most suitable for:

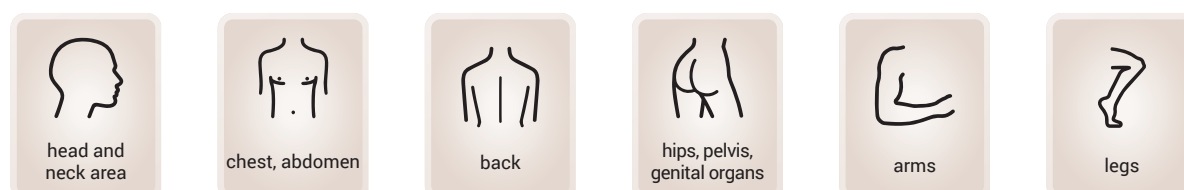

ANELLI set description and use:

ANELLI set is an ideal choice for travellers who always want to have a magnetic therapy device at hand. ANELLI set combines the assets of a small circular CAVO applicator and a small flat MANIC applicator.

The size of MANIC applicator makes it possible for it to generate a strong magnetic field and the applicator is therefore an excellent choice for pain relief. It has also proved to be highly efficient at treating patients with carpal tunnel syndrome, a condition that brings, typically at night, sharp and persistent pain in the wrist and pins and needles in the palms. Pulsed magnetic therapy, if used on a regular basis, is highly indicated in patients with this diagnosis, having the potential to fully substitute standard corticosteroid therapy as well as prevent the necessity to carry out an invasive operation on such a patient. A small MANIC applicator may also be used to treat so called tennis elbows or “frozen shoulder” syndrome. Mothers can use it to treat children with a frequent stomach ache (“baby colic”) or those children suffering from pain as a result of growing teeth.

In contrast, a CAVO applicator creates a homogeneous magnetic field that penetrates the whole depth of the tissue being treated. That is why it is recommended, among others, for alleviating headaches and pain in cervical spine as well as to regenerate joints stricken with arthrosis and rheumatoid arthritis.

ANELLI set applicators:



A circular CAVO applicator

- diameter: 260 mm

Further information about CAVO applicator is available on page 34



A flat MANIC applicator

- dimensions: 170 x 210 mm

Further information about MANIC applicator is available on page 22


ABOLITO SET



ABOLITO set includes:

- a flat applicator PIANO, dimensions: 460 x 490 mm
- a circular applicator CAVO, diameter: 260 mm
- a two-output pulse generator (DUO FORTE)
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas ANELLI set is the most suitable for:

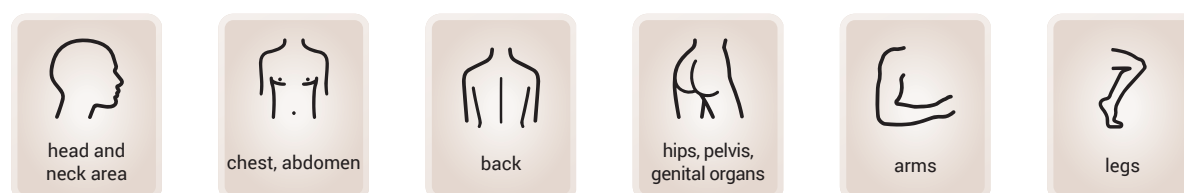

ABOLITO set description and use:

ABOLITO is a magnetic therapy set that includes a circular applicator CAVO and a large flat applicator PIANO.

CAVO is a circular applicator creating a homogenous magnetic field. It is therefore highly suitable for treating tissue structures located deep inside the human body. It is efficient, among other conditions, in treating osteoporosis in upper and lower limbs as well as arthrosis and rheumatoid arthritis. CAVO also substantially accelerates fracture healing processes. As far as arthrosis and rheumatoid arthritis are concerned, CAVO ensures an even distribution of magnetic field to all tissue layers which is a necessary prerequisite for successful treatment of these joint diseases. This makes it possible for ligament and cartilage cells to regenerate, optimizes synovial fluid production and removes swelling and inflammation-related complications. CAVO also provides migraine relief, treats upper respiratory tract, tinnitus and cervical spine pain regardless of origin.

In contrast, a PIANO applicator is indicated in patients with pain in the spinal area and backache in general as well as patients suffering from lumbago and displaced or herniated discs. It has also proved highly efficient in regenerating organs placed in the abdominal and thoracic cavities, in treating problems related to a reduction of blood flow into limbs (especially in patients with diabetes), and in treating inflamed varicose veins.

ABOLITO set applicators:



A flat PIANO applicator

- dimensions: 460 x 490 mm

Further information about PIANO applicator is available on page 25



A circular CAVO applicator

- diameter: 260 mm

Further information about CAVO applicator is available on page 34


MELIO SET



MELIO set includes:

- a flat applicator PIANO, dimensions: 460 x 490 mm
- a flat applicator MANIC, dimensions: 170 x 210 mm
- a two-output pulse generator (DUO FORTE)
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas ANELLI set is the most suitable for:

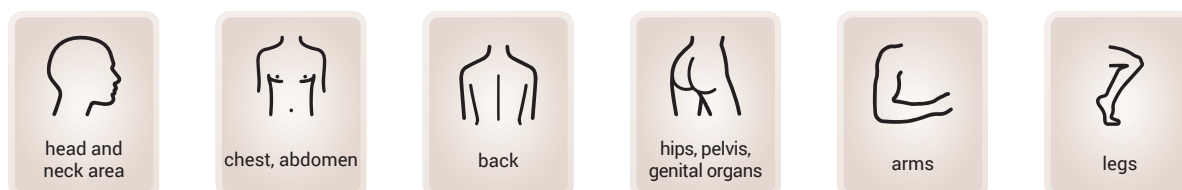

MELIO set description and use:

MELIO is a magnetic therapy set that includes a flat back applicator PIANO and a small flat applicator MANIC.

PIANO is an applicator that may be used, among other conditions, to treat backache as well as pain in any part of spine – cervical, thoracic and lumbar (small of the back). The dimensions of the applicator make it possible for its therapeutic effect to also reach the muscles nearby, whose stiffness and local contractures frequently contribute to the development of the painful condition. PIANO increases blood flow in the exposed tissues, relaxing the muscle fibres and expediting the disposal of acidulous metabolic wastes. It is highly efficient in the chest area as well so it can be used to treat angina pectoris and to regenerate heart muscles after a heart attack. The applicator has also proved to be of great help in treating patients suffering from ulcers in stomach and duodenum, inflammations in pancreas and intestines, constipation, painful menstruation and gynaecological disorders in general.

In contrast, a MANIC applicator is highly suitable for targeted application to places where intensive therapeutic effects are required. It can be used to alleviate pain in joints and teeth, treat ringing in the ear and has proved highly efficient in treating tennis elbow and carpal tunnel syndromes.

MELIO set applicators:



A flat PIANO applicator

- dimensions: 460 x 490 mm

Further information about PIANO applicator is available on page 25



A flat MANIC applicator

- dimensions: 170 x 210 mm

Further information about MANIC applicator is available on page 22






MEDICO SET



MEDICO set includes:

- a flat **applicator PIANO**, dimensions: 460 x 490 mm
- a circular **applicator CAVO**, diameter: 260 mm
- a flat **applicator MANIC**, dimensions: 170 x 210 mm
- a **two-output pulse generator (DUO FORTE)**
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas MEDICO set is the most suitable for:

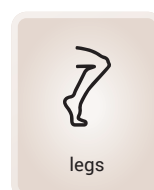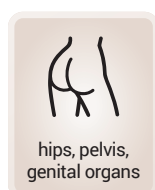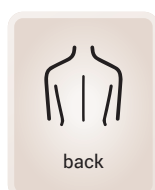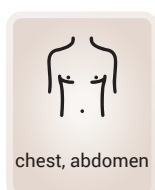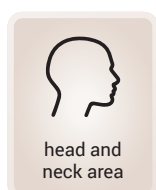



MEDICO set description and use:

MEDICO is a magnetic therapy set that includes a circular applicator CAVO, a wide back applicator PIANO and a small flat applicator MANIC.

PIANO is an applicator that may be used, among others, to regenerate and increase blood flow into lower limbs but is also highly appreciated by people suffering from varicose veins, osteoporosis as well as vascular and neurological complications caused by diabetes. Application to the back helps to alleviate pain in the spinal area as well as to prevent and treat asthma attacks. It is recommended that patients suffering from asthma attacks also apply the device to the chest.

Efficient regeneration of the respiratory tract can also be achieved by using a circular applicator CAVO. This applicator creates a long-range therapeutic magnetic field and, when applied to cervical spine, it is able to treat respiratory tract issues as well. In addition, this is the most suitable applicator to treat all joint diseases.

Finally, MANIC is a small applicator for efficient pain relief and targeted application to e.g. the perineal area (prostate issues, haemorrhoids) and organs placed in the abdominal and thoracic cavities.

MEDICO set applicators:



A flat PIANO applicator

- dimensions: 460 x 490 mm

Further information about PIANO applicator is available on page 25



A circular CAVO applicator

- diameter: 260 mm

Further information about MANIC applicator is available on page 34



A flat MANIC applicator

- dimensions: 170 x 210 mm

Further information about PIANO applicator is available on page 22


LETTINO SET



LETTINO set includes:

- a four-piece **applicator LETTINO**, dimensions: 310 x 1,010 mm
- a flat **applicator MANIC**, dimensions: 170 x 210 mm
- a two-output pulse generator (**DUO FORTE**)
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas LETTINO set is the most suitable for:

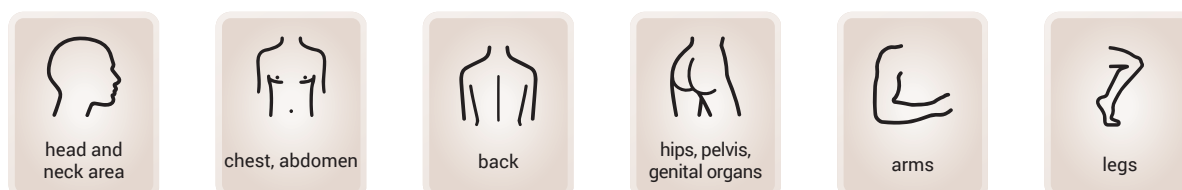

LETTINO set description and use:

This is a magnetic therapy set that includes a medium-sized “bed-type” applicator LETTINO and a small flat applicator MANIC.

The size of a LETTINO applicator covers the whole spinal area including SI joints and pelvis. This is very important because spine works as a whole. Paradoxical as it may seem, a dysfunctional SI joint may cause headaches or cervical spine issues. LETTINO regenerates bones, ligaments and nerves within the whole spinal range, alleviating pain and contributing to restore its proper functions. It can also be used to treat haemorrhoids and prostate issues as well as to fight inflamed varicose veins and deep vein thrombosis. Finally, LETTINO has proved to be efficient in improving blood circulation in limbs.

In contrast, a MANIC applicator has been designed to treat smaller sections of head, abdomen and limbs. When applied to neck, one of the conditions it can be used to treat is hypothyroidism (underactive thyroid). Facial application is an efficient way to treat paradontosis, toothaches and hearing impairments. The applicator is suitable for alleviation of pain in joints and may be safely and efficiently used by children patients as well.

LETTINO set applicators:



A flat LETTINO applicator

- dimensions: 310 x 1,010 mm

Further information about LETTINO applicator is available on page 28



A flat MANIC applicator

- dimensions: 170 x 210 mm

Further information about MANIC applicator is available on page 22






VIVA SET



VIVA set includes:

- a four-piece **applicator LETTINO**, dimensions: 310 x 1,010 mm
- a circular **applicator CAVO**, diameter: 260 mm
- a **two-output pulse generator (DUO FORTE)**
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas VIVA set is the most suitable for:

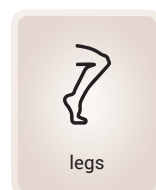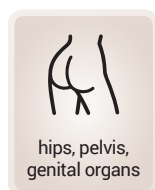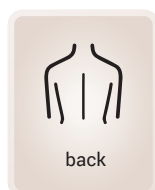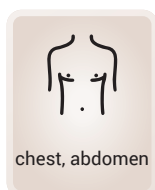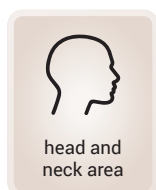


VIVA set description and use:

This is a magnetic therapy set that includes a medium-sized “bed-type” applicator LETTINO and a circular applicator CAVO.

The size of a LETTINO applicator covers the whole spinal area including SI joints and pelvis. This is very important because spine works as a whole. Paradoxical as it may seem, a dysfunctional SI joint may cause headaches or cervical spine issues. LETTINO regenerates bones, ligaments and nerves within the whole spinal range, alleviating pain and contributing to restore its proper functions. It can also be used to treat haemorrhoids and prostate issues as well as to fight inflamed varicose veins and deep vein thrombosis. Finally, LETTINO has proved to be efficient in improving blood circulation in limbs.

CAVO is a circular applicator creating a homogenous magnetic field. It is therefore highly suitable for treating tissue structures located deep inside the human body. It is efficient, among other conditions, in treating osteoporosis in upper and lower limbs as well as arthrosis and rheumatoid arthritis. CAVO also substantially accelerates fracture healing processes. As far as arthrosis and rheumatoid arthritis are concerned, CAVO ensures an even distribution of magnetic field to all tissue layers which is a necessary prerequisite for successful treatment of these joint diseases. This makes it possible for ligament and cartilage cells to regenerate, optimizes synovial fluid production and removes swelling and inflammation-related complications. CAVO also provides migraine relief, treats upper respiratory tract, tinnitus and cervical spine pain regardless of origin.

VIVA set applicators:

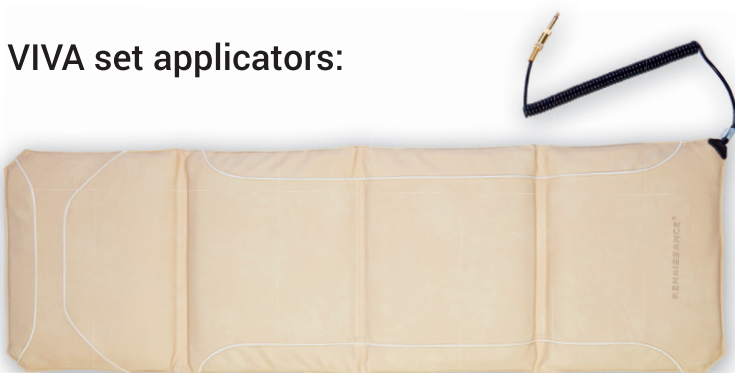

A flat LETTINO applicator

- dimensions: 310 x 1,010 mm

Further information about LETTINO applicator is available on page 28



A circular CAVO applicator

- diameter: 260 mm

Further information about CAVO applicator is available on page 34






VIVA PLUS SET



VIVA PLUS set includes:

- a four-piece **applicator LETTINO**, dimensions: 310 x 1,010 mm
- a circular **applicator CAVO**, diameter: 260 mm
- a flat **applicator MANIC**, dimensions: 170 x 210 mm
- a **two-output pulse generator (DUO FORTE)**
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas VIVA PLUS set is the most suitable for:

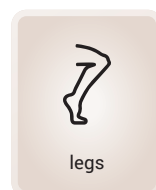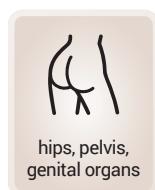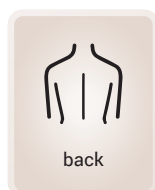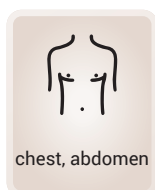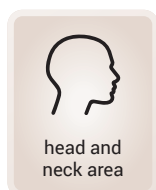

VIVA PLUS set description and use:


This is a magnetic therapy set that includes three applicators: a medium sized “bed-type” applicator LETTINO, a circular applicator CAVO and a small flat applicator MANIC. This combination makes it possible for you to make the most of what low-frequency pulsed magnetic therapy offers.

The size of a LETTINO applicator covers the whole spinal area including SI joints and pelvis. This is very important because spine works as a whole. Paradoxical as it may seem, a dysfunctional SI joint may cause headaches or cervical spine issues. LETTINO regenerates bones, ligaments and nerves within the whole spinal range, alleviating pain and contributing to restore its proper functions. It can also be used to treat haemorrhoids and prostate issues as well as to fight inflamed varicose veins and deep vein thrombosis. Finally, LETTINO has proved to be efficient in improving blood circulation in limbs.

CAVO is a circular applicator creating a homogenous magnetic field. It is therefore highly suitable for treating tissue structures located deep inside the human body. It is efficient, among other conditions, in treating osteoporosis in upper and lower limbs as well as arthrosis and rheumatoid arthritis. CAVO also substantially accelerates fracture healing processes. As far as arthrosis and rheumatoid arthritis are concerned, CAVO ensures an even distribution of magnetic field to all tissue layers which is a necessary prerequisite for successful treatment of these joint diseases. This makes it possible for ligament and cartilage cells to regenerate, optimizes synovial fluid production and removes swelling and inflammation-related complications. CAVO also provides migraine relief, treats upper respiratory tract, tinnitus and cervical spine pain regardless of origin.

The size of the MANIC applicator makes it possible for it to generate an intensive magnetic field which can alleviate a number of kinds of pain. Its analgesic and regenerating effect is highly appreciated by people suffering from arthrosis or rheumatoid arthritis in knee and hip joints as well as patients with tennis elbows and carpal tunnel syndrome. It can also be used for intensive highly targeted therapy and regeneration of various organs and tissues as well as to get patients rid of muscular seizures. And it is of equally great help in patients suffering from ringing in the ear (tinnitus, Ménière's disease) and otitis media. Finally, the applicator is also highly recommended for patients with toothache and paradontosis.

VIVA PLUS set applicators:

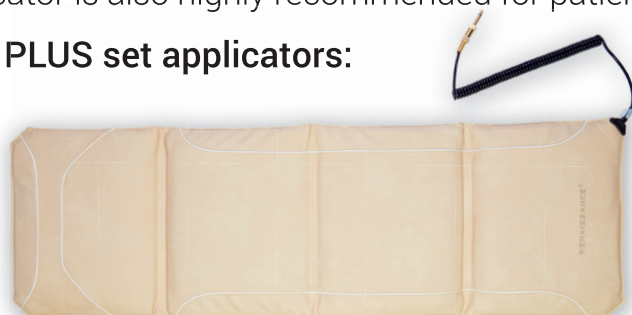

A flat LETTINO applicator

- dimensions: 310 x 1,010 mm

Further information about LETTINO applicator is available on page 28

VIVA PLUS set applicators:



A flat MANIC applicator

- dimensions:
170 x 210 mm

Further information about MANIC applicator is available on page 22



A circular CAVO applicator

- diameter: 260 mm

Further information about CAVO applicator is available on page 34


SIECCO SET



SIECCO set includes:

- a five-piece **applicator CASSO**, dimensions: 550 x 1,730 mm
- a flat **applicator MANIC**, dimensions: 170 x 210 mm
- a **two-output pulse generator (DUO FORTE)**
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas SIECCO set is the most suitable for:

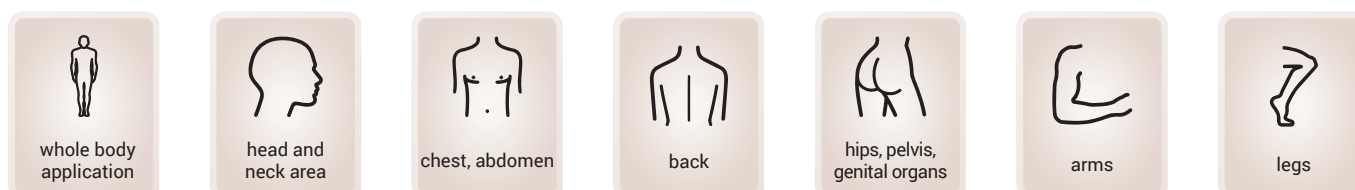

SIECCO set description and use:

This magnetic therapy set includes a large “bed-type” applicator CASSO and a small flat applicator MANIC.

The size of a CASSO applicator makes it an ideal choice for whole body regeneration. Lie down and enjoy the analgesic, anti-inflammatory, and blood circulation enhancing effects of this low-frequency pulsed magnetic therapy device from head to toe. CASSO is designed in such a way that ensures magnetic field intensity is equally distributed to the whole area of the applicator. This means each part of your body is provided sufficient treatment. This is incredibly important because human musculoskeletal system works as an inseparable whole. Headaches, to provide just one example, are often caused by dysfunctional SI joints, gluteal muscles or even knees. If this is the case, the CASSO is a great choice as it is able to regenerate each of the body sections, addressing the actual cause of the problems.

In contrast, the shape of a MANIC applicator makes it possible for it to generate strong highly-focused magnetic field. As a result, it is highly suitable for targeted regeneration of various body parts as well as to alleviate pain. It can be used to treat a number of injury types as well as toothache. And it is also recommended to alleviate pain in joints and spinal sections.

SIECCO set applicators:



A flat CASSO applicator

- dimensions: 550 x 1,730 mm

Further information about CASSO applicator is available on page 31



A flat MANIC applicator

- dimensions: 170 x 210 mm

Further information about MANIC applicator is available on page 22


VIERA SET



VIERA set includes:

- a circular applicator VIERA, diameter: 650 mm
- a two-output pulse generator (DUO FORTE)
- power supply
- a magnetic field tester
- a user's manual, information and picture guide
- a quality case to safely store and transport the device

Body areas VIERA set is the most suitable for.

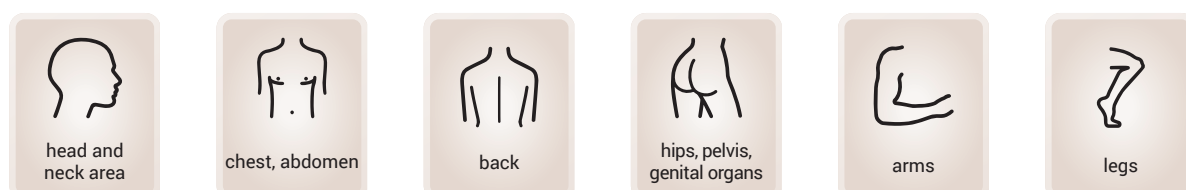

VIERA set description and use:

VIERA is a magnetic therapy set that includes an eponymous applicator. The set is primarily designed for professional spa and physical medicine/rehabilitation facilities as well as physiotherapists' offices and, last but not the least, healthcare providing homes and day-care centres for elderly and disabled people.

The VIERA applicator makes it possible for a patient to be placed inside a homogeneous (direct-current) pulsed magnetic field that evenly penetrates the entire depth of the tissues being treated. The applicator's regenerating, analgesic and anti-inflammatory effects are fully indicated in patients suffering from arthrosis in hip and knee joints, rheumatoid arthritis, backache regardless of origin, osteoporosis and a number of other diseases. While pulsed magnetic therapy by means of this applicator is an efficient treatment method to fight all the foregoing diseases, it is especially recommended to treat osteoporosis, a disease whose riskiest symptom is decreased bone density in the femoral neck area and in the area of lumbar vertebrae.

VIERA set applicator:



A circular VIERA applicator

- diameter: 650 mm

Further information about VIERA applicator is available on page 39



WHEN MAGNETIC THERAPY SHOULD NOT BE USED OR CAUTION EXERCISED

When care and caution should be exercised


- Patients with a cochlear implant may not apply magnetic therapy to head, neck or upper torso
- Patients with hyperfunction (increased function) of endocrine glands do not apply magnetic therapy devices to the respective area.
- Patients diagnosed with a tumour, patients who have had a cancer surgery, and patients on chemotherapy or radiotherapy should consult their physician prior to application.
- An infectious disease during which body temperature exceeds 38°C (100.4°F).
- A doctor's consent is required in patients with active tuberculosis.
- Patients with neurological disorders accompanied by seizures (e.g. epilepsy) should not apply magnetic therapy devices to the head area (unless a neurologist recommends otherwise).
- In patients with serious mycotic (fungal) diseases.
- Combining magnetic therapy and radiotherapy is not recommended.

When it should not be used


- In patients who use sensitive electronic equipment such as cardiac pacemakers, implantable cardioverter defibrillators, etc.
- In pregnant women. Although there is no evidence magnetic field negatively affects pregnant women or unborn children, caution is highly advised.
- In patients suffering from haemorrhage. Magnetic field increases blood circulation in tissues which may, in some cases, exacerbate the condition.
- If you have any questions related to pulsed magnetic therapy treatment, please, consult your physician!

MAGNETIC THERAPY APPLICATION POINTS

The following human body drawings include numbers that indicate areas to which magnetic field applicators should be applied in order to treat particular diseases (see pp. 67 and 68). These application points are the result of long-term research carried out by doctors from both the Czech Republic and abroad.

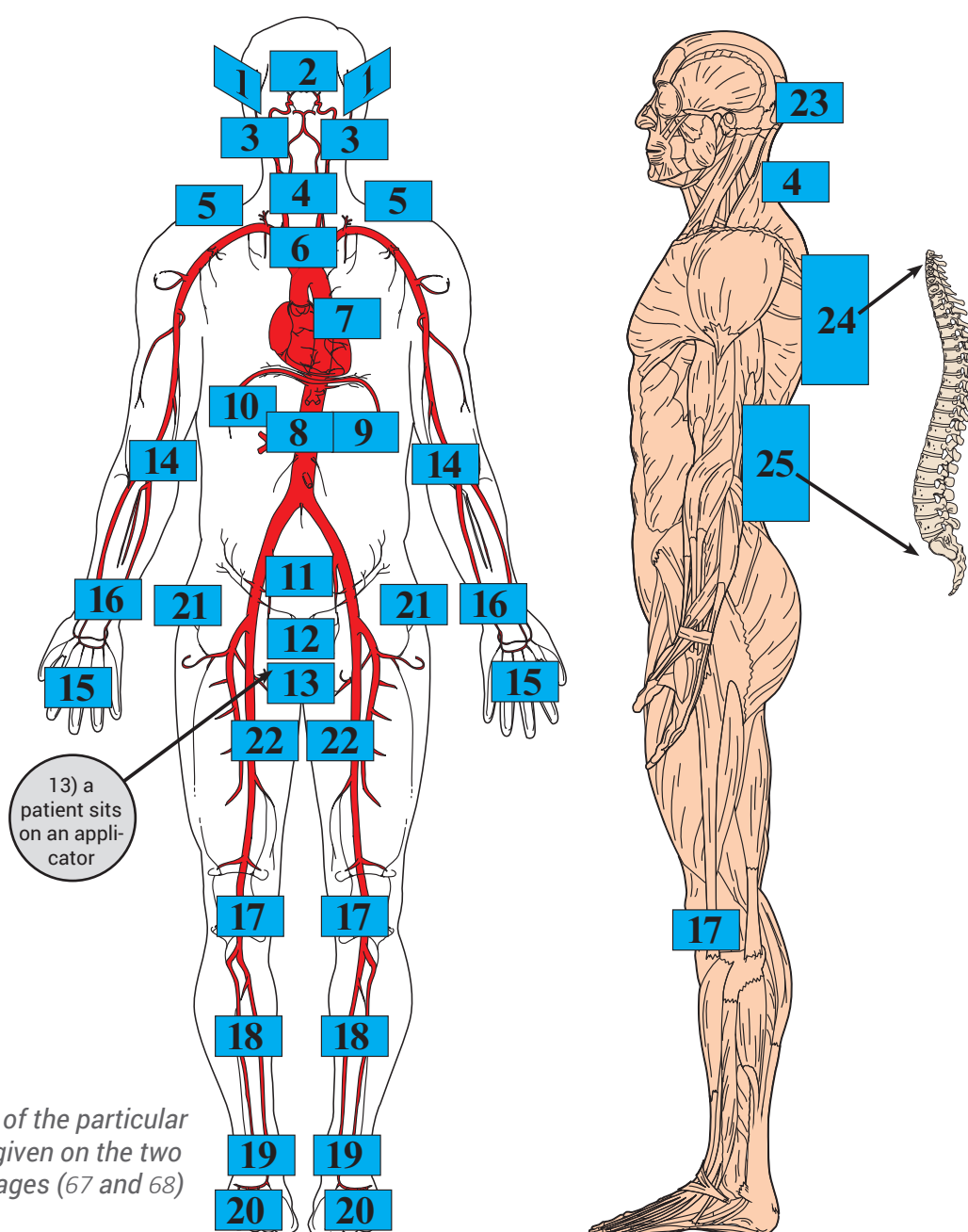

Description of the particular points are given on the two following pages (67 and 68)



1) otitis media (in the right, left or both ears), initial phase when exudate caused by inflammation is missing

mandibular (jaw) joint pain

tinnitus (ringing in the ear),
medication supported

Parkinson's syndrome,
medication supported

migraines

2) inflammation in frontal
sinuses (initial phase)

common and chronic cold

retinal regeneration

optic atrophy, as supplementary
medication

corneal perforation (injuries, surgeries),
helps during final healing stages

3) pain and inflammation in teeth

paradontosis

post-extraction pain (pain experienced
after tooth extraction, though not
immediately after extraction!)

trigeminal nerve pain, as
supplementary medication

4) pain caused by blocked cervical spine

headaches

carpal tunnel syndrome

impaired blood circulation in
arms (pins and needles)

cold – hoarseness

pharyngitis

5) stiff painful shoulders
(so called frozen shoulders)

arthralgia or arthrosis in shoulder
joints (pain in joints)

6) asthma, as supplementary medication
bronchitis, as supplementary medication

high blood pressure

hiccups

7) heart muscle regeneration once
a heart attack has subsided

after by-pass surgery

coronary artery disease
(ischemic heart disease)

tachycardia (fast heart rhythm)

8) inflammation and spasms in
stomach, sour stomach

chronic gastric or duodenal
ulcer disease (not to be applied
at the time bleeding occurs)


9) inflammation in pancreas –
pancreatitis (chronic)

diabetes (applicator is to be
applied to the left lower side)

10) liver regeneration

11) female sterility caused by
inflamed fallopian tubes

constipation, problems with
bowel movement

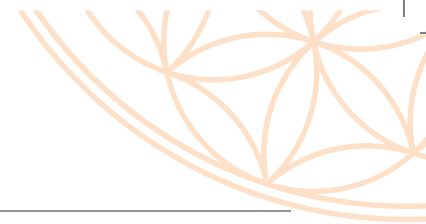

12) cystitis (bladder infection), particularly in women, as supplementary treatment
prostatitis in men
numbers 12 a 25 – urinary incontinence in children (enuresis)

13) prostatitis
male sexual dysfunction
erectile dysfunction
haemorrhoids (unless bleeding is present)

14) acute and chronic “tennis elbow”, elbow pain (arthrosis)

15) inflammation in joints in final healing stages
rheumatic diseases
polyarthritis progresiva (rheumatoid arthritis) in wrist

16) numbers 16 and 4 – carpal tunnel syndrome

17) pain to arthrosis in knee joints
poor blood circulation in lower limbs in patients with diabetes (diabetic leg)

18) poor blood circulation in lower limbs in patients with diabetes (diabetic leg)
osteoporosis in shin-bone
chronic varicose ulcers (unless there is bleeding), as supplementary medication
inflammation of bone marrow (osteomyelitis) together with antibiotics!

19) pain in insteps and ankles
swollen ankles
poor blood circulation in lower limbs in patients with diabetes (diabetic leg)

20) rheumatic diseases
gout, painful heel spurs
poor blood circulation in lower limbs in patients with diabetes (diabetic leg)

21) pain and arthrosis in hip joints
osteoporosis in femoral neck area

22) poor blood circulation in lower limbs in patients with diabetes (diabetic leg)

23) migraines
Parkinson's syndrome, medication supported
spasms (non-epileptic)
phantom pain in missing limbs

24) regeneration of cervical and thoracic spine
carpal tunnel syndrome
partial hearing loss
degeneration of spinal chord of various origins
spinal apophyses (outgrowths)
Bechterew's disease
osteoporosis in spine
phantom pain in missing limbs

25) sciatica, pain in loins and small of the back
painful kidneys, helps to flush out kidney stones at early stages of kidney stone disease
osteoporosis in spine



CONTACT INFORMATION



Nice Care
Health Products Limited

web: www.nicecare.uk


email: info@nicecare.uk

phone: 0203 488 6267

Nice Care Health Products Limited is a UK registered company
Company number: 12667901





NOTES

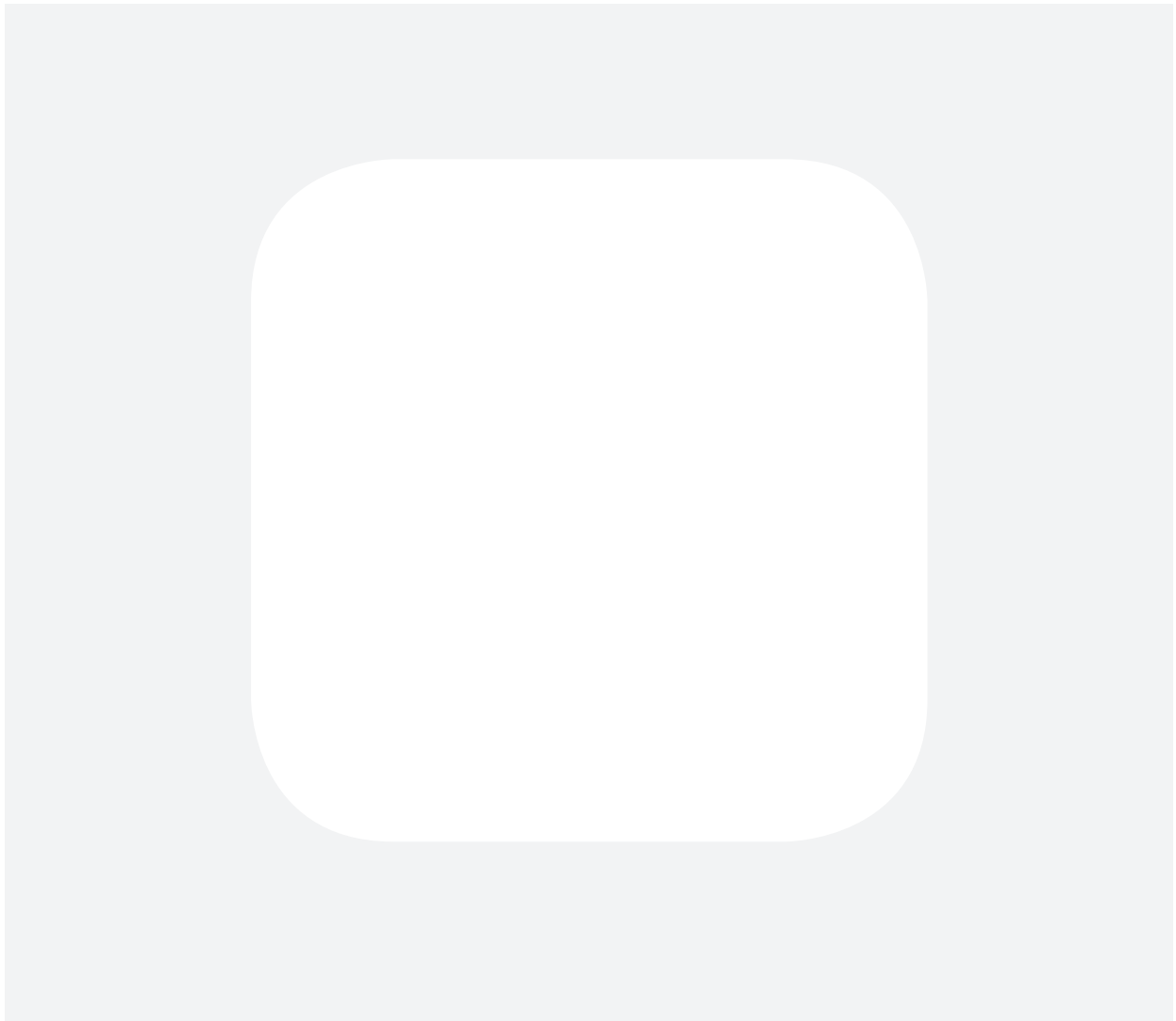



NOTES

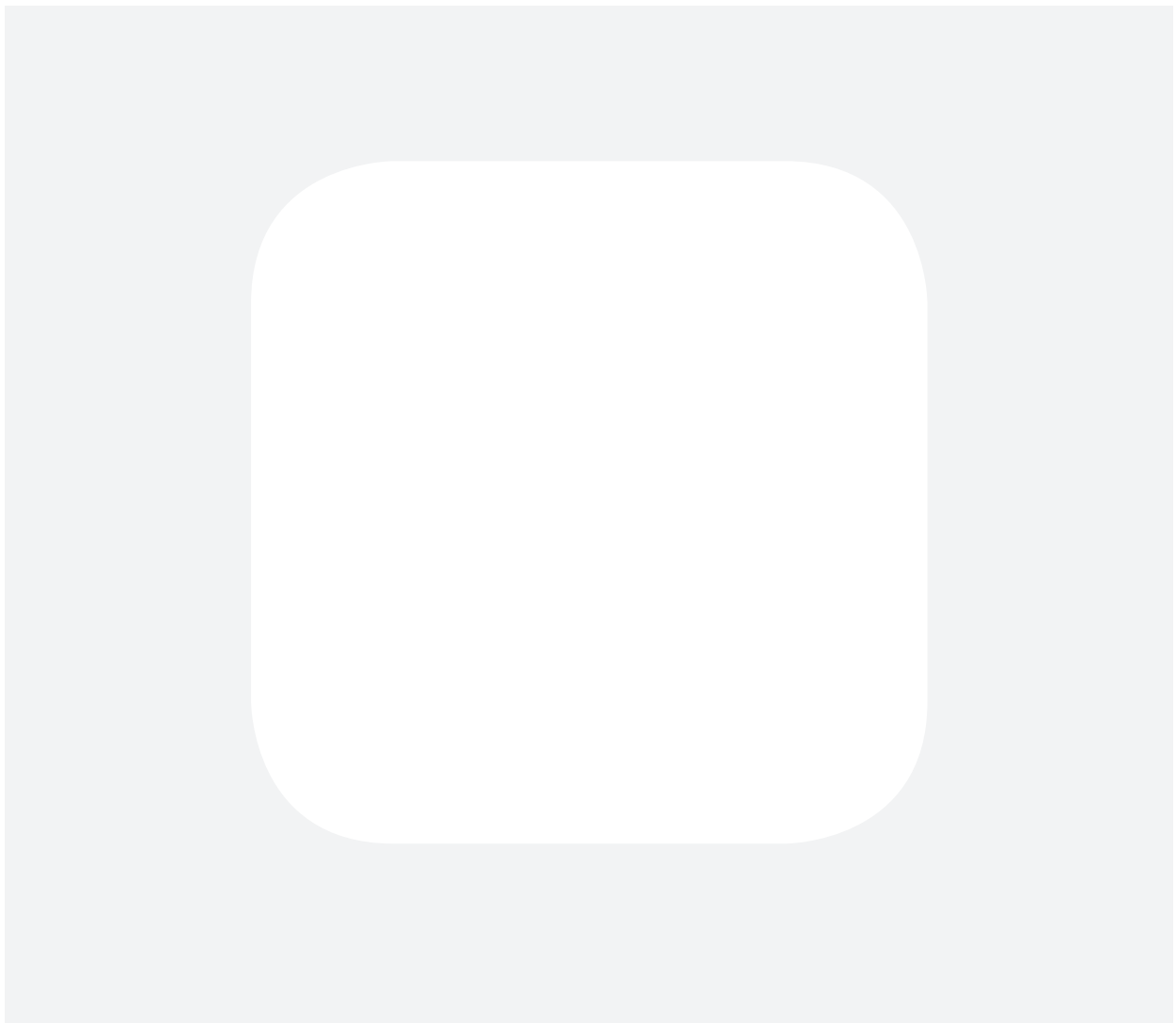

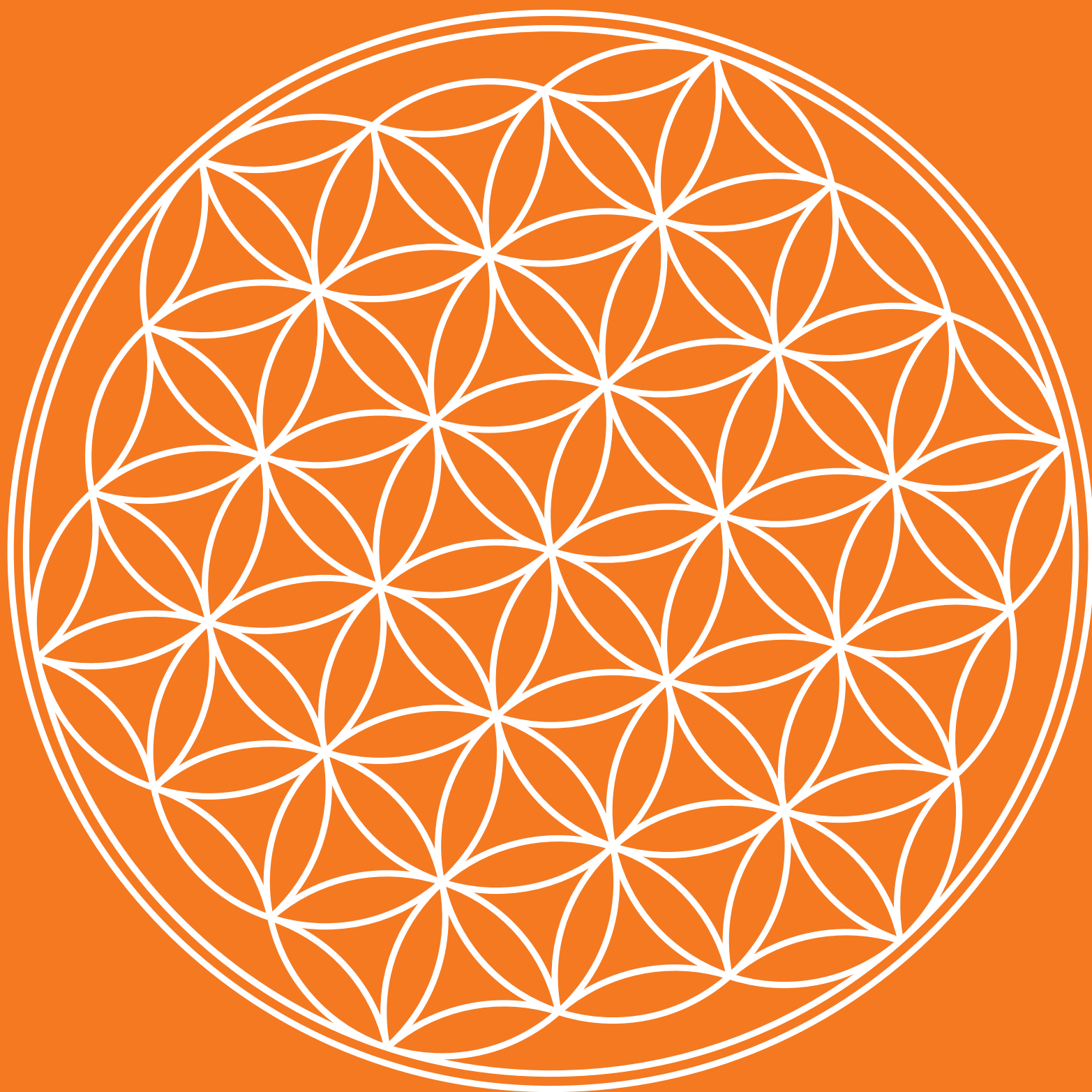

RENAISSANCE[®]

Live a life without pain or limitations



Nice Care
Health Products Limited